



Beaver-Visitec® Dry Eye Products

Discover why BVI is the world leader in punctal occlusion — providing solutions for every stage in your patients' continuum of care.



Permanent



Long-Term



Mid-Term




Short-Term



Partial



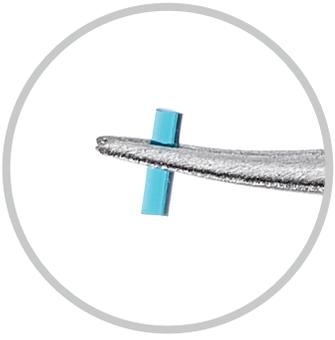

Mid-Term Occlusion Extend® Absorbable Synthetic Implants

Extend absorbable synthetic implants provide occlusion for up to 3 months for patients experiencing short-term dry eye symptoms. These implants are ideal for the transient dry eye associated with LASIK or other ocular surgeries. Extend implants can also be used for allergy sufferers, contact lens discomfort, climate change or to sustain drug contact time and effectiveness.




Short-Term Occlusion Collagen Implants

These comfortable and economical implants provide occlusion for 7-10 days. They are ideal for determining if permanent punctal occlusion is right for your patient.



Partial Occlusion Micro Flow™

Designed to provide a simple, reversible means of partial occlusion, micro flow punctal occluders limit the drainage of tears into the lacrimal duct. Ideal when total occlusion is too much.

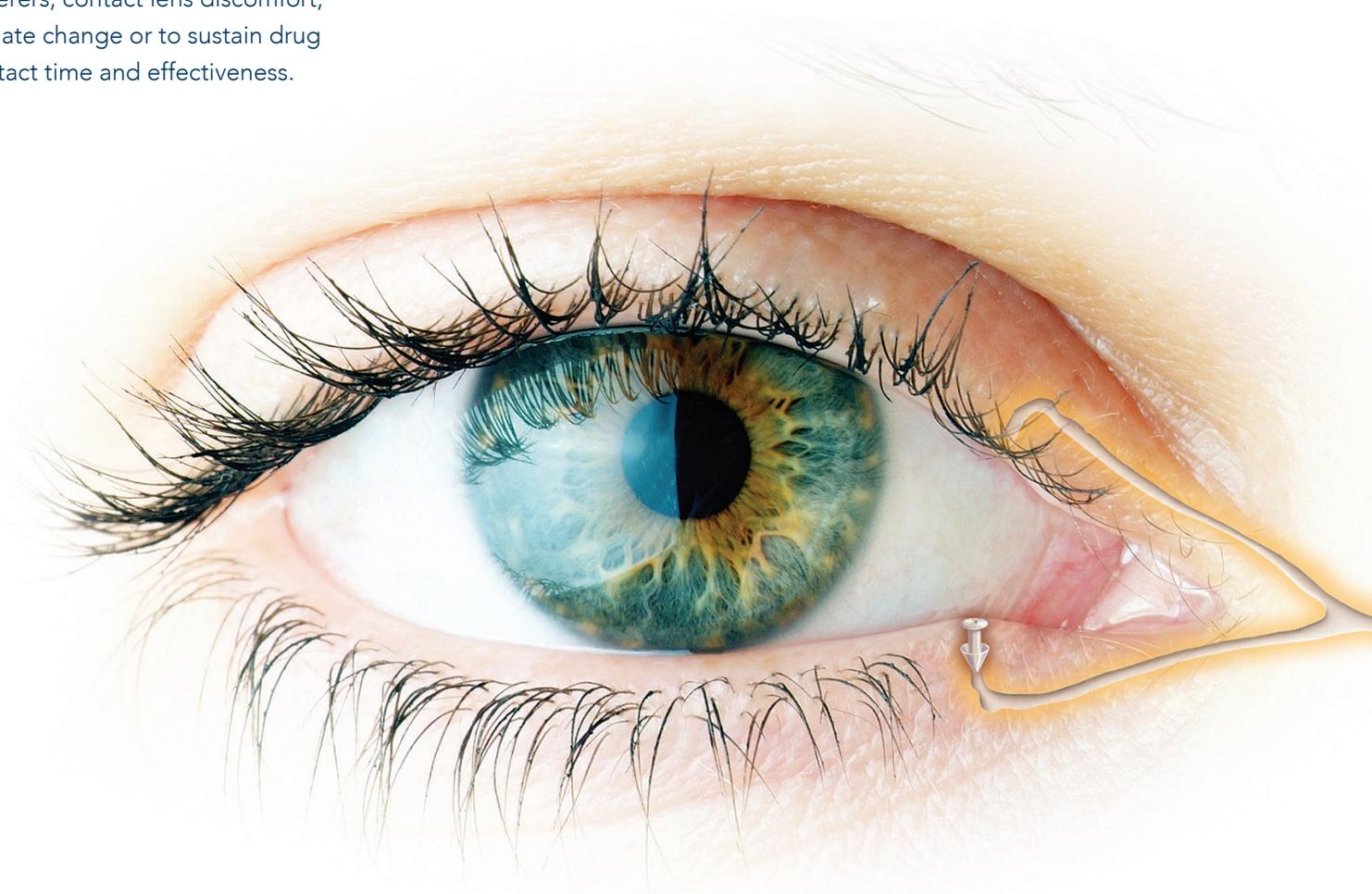



Permanent Occlusion Parasol® Punctal Occluder

Permanent occlusion, though easily reversible, is indicated for chronic dry eye conditions. Made of medical grade silicone, one glance at parasol punctal occluder and you'll recognize a breakthrough product, offering distinct benefits for you and your dry eye patients.

Easy To Insert


The hollow nose collapses on insertion, eliminating the need for punctal dilation in most cases. It's easier on you and more comfortable for your patients.

Simple Sizing

Two sizes, small and medium, fit virtually all patients. X-small and large sizes are available for unusually sized puncta.

Guaranteed Retention*

Parasol is designed for excellent retention, as demonstrated by an independent clinical study*. Even better, parasol is backed by a no-nonsense retention guarantee. If your patient loses our plug within 30 days of insertion, we will replace it at no charge.

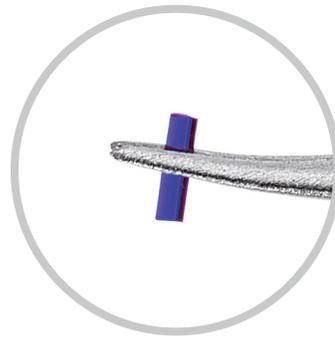

Long-Term Occlusion Extend 180® Absorbable Synthetic Implant

Extend 180® absorbable synthetic implants provide 6 months of effective dry eye relief for patients. These implants can be used for post-surgical or seasonal dry eye relief, contact lens intolerance or dry eye associated with digital eye strain.



Parasol Size Comparison Chart

Traditional Plugs	0.2 mm	0.3 mm	0.4 mm	0.5 mm	0.6 mm	0.7 mm	0.8 mm	0.9 mm +
Parasol®	X-SMALL	SMALL		MEDIUM	LARGE			

Small and medium sizes should accommodate 95% of all patients.
BVI recommends undersizing when in doubt.

* "Plugs Reduce Dry Eye Symptoms, Improve Vision", Craig McCabe, M.D., Review of Ophthalmology, 2009.

Permanent Occlusion

Parasol® Sterile Pre-Loaded Punctal Occluders		Item #	Quantity
	Occluders are packaged in a convenient (1) pair tray. Each occluder is individually sterile and pre-loaded on a high-quality disposable inserter.		
	X-Small (0.25 mm — 0.35 mm)	A14-200	1 pair tray; 2 occluders total
	Small (0.35 mm — 0.65 mm)	A14-201	
	Medium (0.6 mm — 0.85 mm)	A14-202	
Large (0.9 mm+)	A14-203		
Parasol® Sterile Quickload™ Punctal Occluders		Item #	Quantity
	Combines sterile convenience with bulk economy. Occluders are supplied sterile, pre-mounted on a stand. Reusable inserters sold separately.		
	Small (0.35 mm — 0.65 mm)	A14-201QL	8 pair per box; 16 occluders total
	Medium (0.6 mm — 0.85 mm)	A14-202QL	
	Large (0.9 mm+)	A14-203QL	
Combo Pack (4 pair small, 4 pair medium)	A14-205QL		
Parasol® Non-Sterile Bulk Punctal Occluders		Item #	Quantity
	Occluders are supplied non-sterile; 2 plugs per poly bag.		
	Small (0.35 mm — 0.65 mm)	A14-201BN	2 plugs per poly bag; 10 pair per box; 20 occluders total
	Medium (0.6 mm — 0.85 mm)	A14-202BN	
	Large (0.9 mm+)	A14-203BN	
Combo Pack (5 pair small, 5 pair medium)	A14-209BN		

Long-Term Occlusion

Extend 180® Synthetic Absorbable Implants		Item #	Quantity
	0.3 mm X 2.0 mm single pair	A13-2033SP	1 pair box; 2 occluders total
	0.4 mm X 2.0 mm single pair	A13-2044SP	
	0.5 mm X 2.0 mm single pair	A13-2055SP	
	0.3 mm X 2.0 mm bulk pack	A13-2033	10 pair box; 20 occluders total
	0.4 mm X 2.0 mm bulk pack	A13-2044	
	0.5 mm X 2.0 mm bulk pack	A13-2055	
	Extend 180 Absorbable Synthetic Implant variety pack	A13-2000	10 pair box includes: 2 pairs — 0.3 mm X 2.0 mm 5 pairs — 0.4 mm X 2.0 mm 3 pairs — 0.5 mm X 2.0 mm

Mid-Term Occlusion

Extend® Synthetic Absorbable Implants		Item #	Quantity
	0.2 mm x 2.0 mm	A13-202	20 implants per box; packaged in 10 individual pairs
	0.3 mm x 2.0 mm	A13-203	
	0.4 mm x 2.0 mm	A13-204	
	0.5 mm x 2.0 mm	A13-205	
	Absorbable Synthetic Implants (variety pack)	A13-200	10 pair box includes: 3 pair — 0.3 mm 5 pair — 0.4 mm 2 pair — 0.5 mm

Short-Term Occlusion

Collagen Implants		Item #	Quantity
	0.2 mm x 2.0 mm — 10 pair box	A12-102	60 implants per box; packaged in 10 foam sleeves with 6 implants each
	0.3 mm x 2.0 mm — 10 pair box	A12-103	
	0.4 mm x 2.0 mm — 10 pair box	A12-104	
	Variety Pack — 18 implants — 0.2 mm 24 implants — 0.3 mm 18 implants — 0.4 mm	A12-200	

Partial Occlusion

Micro Flow™ Sterile Pre-Loaded Punctal Occluders		Item #	Quantity
	Individually sterile; pre-loaded onto a high-quality disposable insertion instrument		
	Small (0.4 mm — 0.55 mm)	A14-401	1 pair tray; 2 occluders total
	Medium (0.55 mm — 0.7 mm)	A14-402	
	Large (0.7 mm — 0.85 mm)	A14-403	
Micro Flow Non-Sterile Bulk Punctal Occluders		Item #	Quantity
	Occluders are supplied non-sterile; 2 plugs per poly bag.		
	Small (0.4 mm — 0.55 mm)	A20-401BN	5 pair non-sterile bulk pack; 10 occluders total
	Medium (0.55 mm — 0.7 mm)	A20-402BN	
	Large (0.7 mm — 0.85 mm)	A20-403BN	

Other Dry Eye Products


Other Dry Eye Products		Item #	Quantity
Reusable Insertion Instrument Non-sterile reusable insertion instrument for use with bulk punctal occluders		A14-300	1 pair
Punctal Gauge/Micro Dilator Dual purpose instrument intended for determining the proper size parasol to utilize and dilating extra small puncta. Use with bulk punctal occluders.		A14-302	1 instrument
Jewelers Forceps #7 Curved tip forceps are ideal for insertion of collagen and extend implants, as well as removal of the parasol, if needed.		3740	1 instrument


Dry Eye Diagnostic Tests		Item #	Quantity
	BioGlo™ Fluorescein	A17-502	Box of 100 strips
		A17-505	Box of 300 strips
	Lissamine™ Green	A17-501	Box of 100 strips
	Zone-Quick™	A17-503	Box of 50 threads
	Tear Flo™	A17-504	Box of 100 strips

Resource Materials

Billing guidelines, wall charts, and patient brochures are available upon request; Instructional video, dry eye questionnaires, and procedure consent forms are also available.



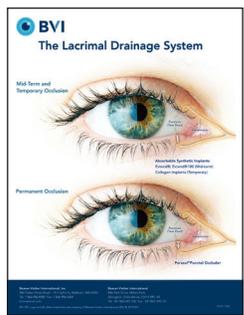

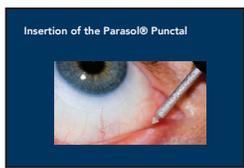

Call your local sales rep or customer service at 1-866-906-8080

Beaver-Visitec International, Inc.
 500 Totten Pond Road – 10 CityPoint, Waltham, MA 02451
 Tel: 1.866.906.8080 Fax: 1.866.906.4304
 bvimedical.com

Beaver-Visitec International Sales Ltd
 85b Park Drive, Milton Park,
 Abingdon, Oxfordshire, OX14 4RY, UK
 Tel: 44.1865.601.256 Fax: 44.1865.595.761

